


# CARONY

## the transport wheelchair

Easy and comfortable – all the way into the car

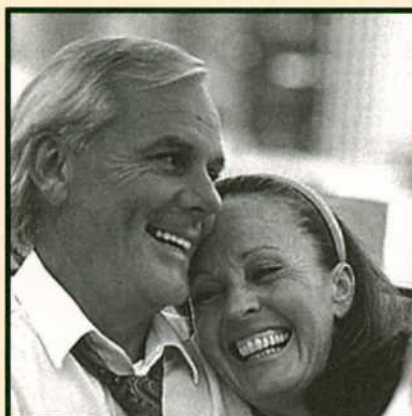

Mobility Conversions Ltd, Langley Road South, Salford, Manchester, M6 6SY

Tel: 0161-745-9737 or Fax: 0161-743-0190



Email: [info@mobilityconversions.co.uk](mailto:info@mobilityconversions.co.uk)



Web Site: [www.mobilityconversions.co.uk](http://www.mobilityconversions.co.uk)



# wheelchair that increases your freedom!

Mobility Conversions Ltd, Langley Road South, Salford, Manchester, M6 6SY



Tel: 0161-745-9737 or Fax: 0161-743-0190



Email: [info@mobilityconversions.co.uk](mailto:info@mobilityconversions.co.uk)

Web Site: [www.mobilityconversions.co.uk](http://www.mobilityconversions.co.uk)

## RIGHT DRIVE POSITION

The Carony chair has an ergonomic tested drive position for better manoeuvrability. The chair is also adjustable in height.

## SUBSTANTIAL DRIVE GRIP

The Drive Grip is ergonomically adjusted and easy to unload and put on.

## ARMREST

The Carony chair has substantial and adjustable armrests as standard.

## VERTICALLY ADJUSTABLE WHEELCHAIR FRAME

The wheelchair frame is vertically variable with 115 mm to be conformed for different car models and ground conditions when you dock into the car. You can vertically adjust the Carony with an easy manoeuvrable crank on the side of the wheelchair frame.



## SECURITY FIRST

The chair has side support and a headrest. It is also crash-tested with the actual chair in the car and meets the security requirements from all the major car manufacturers. Carony is equipped with security functions to minimise risks when loading and unloading. Carony Transport Wheelchair is CE marked.

## COMFORTABLE CHAIR

The anatomically constructed chair has good side support, sway-back support and correct sit angle.

## FREE FOOT PASSAGE

The U-shaped wheelchair frame gives free passage for the passenger's feet when the frame is removed.

## SUBSTANTIAL WHEELS

The substantial wheels are compact and suitable for all environments.

## Simple to manoeuvre, drive and load.



4



5



6

Release the safety stop on the wheelchair frame and push the chair into the car. Unload the wheelchair frame from the swivel seat in the car.

Turn the chair into position, be sure that the turn stop locks and fasten the seat belt. The transfer is now complete!

Store the wheelchair in the trunk and make sure it is well anchored.

# CARONY - the transport

Mobility Conversions Ltd, Langley Road South, Salford, Manchester, M6 6SY



Tel: 0161-745-9737 or Fax: 0161-743-0190



Email: [info@mobilityconversions.co.uk](mailto:info@mobilityconversions.co.uk)

Web Site: [www.mobilityconversions.co.uk](http://www.mobilityconversions.co.uk)

Carony is a unique, manageable transport wheelchair system that eliminates lifting people and simplifies transferring into and out of private vehicles. With assistance from Carony, you can independently move your passenger without straining your arms or back!

The Carony system consists of a wheelchair frame, a comfortable seat that you use both in and outside of the vehicle and a mounted swivel seat base. With a couple of effortless manoeuvres, you are able to position the chair and passenger to the front passenger space. The wheelchair frame can then be stored in the trunk as it takes up minimal space.

Once in place the passenger will use the original seat belt. The chair has been especially designed, crash tested and approved, meeting all valid legislation. Security and comfort is equal or better to the original seat and practical, optional accessories make the chair flexible to suit to any special needs or demands. Carony also has the flexibility to load when the passenger chooses to travel by public transport.

Carony will fit most cars and MPVs regardless of size or number of doors. Approved AutoAdapt-BEV dealers will have fully trained personal with responsibility for function, security and installation of the Carony Transport Wheelchair System.



## Comfortable transportation without lifting!



*Carony is so simple to drive and manoeuvre, especially in small areas. The ergonomic shape gives an excellent driving position.*



*The chair is easy to transport from both inside and out of the car; you start by turning the swivel seat out.*



*Remove the drive grip and footrest and dock the chair to the swivel seat in the car.*